


# Elongated swab

X-RAY

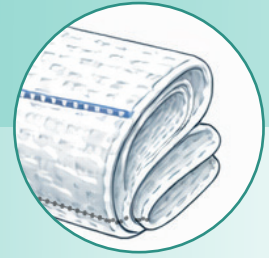

## A unique innovation on the market!

- ✓ Solution designed for minimally invasive surgery
- ✓ Compatible with **laparoscopic** trocars
- ✓ **Elongated** shape for easier insertion



REF 070315

**15 cm**

Trocars 5 mm  
(width: 1 cm)

REF 070320

**20 cm**

Trocars 10 mm  
(width: 2.5 cm)

*Handle the swab by its end during insertion and removal through the trocar.*

## Technical specifications

- **50 g non-woven material:** high absorption capacity
- **X-RAY detectable:** enhanced radiographic identification for precise visualization
- **Swabs stitched along the length:** increased strength and secure fixation of the barium thread
- **Innovative design:** easy passage through a 5 mm trocar
- **Available in two sizes: 15 cm (slim) for 5 mm trocars and 20 cm (thicker) for 10 mm trocars,** depending on your absorption requirements
- **Packs of 5 sterile swabs in double sterile packaging:** optimal safety in the operating room
- **Packaging:** carton containing 4 boxes of 50 pouches

Cécile LECLERC-CHALVET, Product Manager  
c.leclerc@raffin-medical.com +33 (0)4 27 29 00 15

Designs and models published in the French Official Bulletin of Industrial Property (BOPI) No. 26/05 of April 24, 2026.  
This is a Class IIa medical device manufactured by our subsidiary V-MEDICAL PRODUCTS. This medical device is a regulated healthcare product and bears the CE0197 marking in accordance with applicable regulations. Its European Authorized Representative is RAFFIN MEDICAL.

